



BABY SCALE **BAK**
User Manual



v.201811

Value Each Gram

1. Technical specifications

- Work temperature: 10 ~ 40°C
- Environmental humidity: <80%
- Maximum capacity: 20kg
- Graduation: 5g
- Power supply: AA X 4 batteries
- Low power: L o
- Error of indication: Err
- Unit change: kg or lb
- Memory function: optionally save and recall weighing value
- Turn off automatically

2. Keys' Functions

- **ON/ZERO(/OFF):**
Press the switch to power on the scale and start up weighing process. For zero operation, press the switch again to clear the display and "zero" annunciator will show in the left low corner of LCD. To turn off the scale please hold pressing this switch for about 3 seconds.
- **MODE**
Weighing units exchange switch. Press the switch to select kg and lb in turn. The suitable indicator of unit is shown in LCD right side by kg or lb annunciator.
- **MEM:**
While the scale is in "0.00" situation, press this switch to recall the former one weighing value saved in the memory and the suitable indicator is lighted. During weighing if the switch is pressed and the stable value will be saved, otherwise twinkling value is not saved in the memory.


3. Operation

- After the batteries placed, press "**ON/ZERO (/OFF)**" key to start the scale until "**0.00**" appearance.
- Gently place baby on the pan, weighing record will be flashed on the window after a few seconds. If baby remain lying on the pan after its weight value has shown, the weighing value will continue to flash for 5 minutes and then turn off automatically.
- If take baby from the pan immediately after baby's weight has been flashed, the scale will return to "**0.00**" and then turn off automatically within 5 minutes if no operation again.
- If see about the former saved weighing value, please press the "**MEM**" switch while the scale is in "**0.00**" position. The scale will save optional value during each weighing operation.

4. Cautions

- *The scale should be placed on a flat and steady surface for baby's easement and safety*
- *You must replace batteries immediately when window display "L o"*
- *Keep the scale clean to ensure the normal working.*
- *Take the batteries out if no operation for a long term.*
- *Please place baby in the center of the pan otherwise probably affect weighing accuracy.*
- *Overload (heavier than 23.000kg) is strictly prohibited.*

Value Each Gram

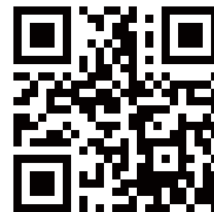

v.201811

User Manual **BAK**